


Balancium

Guide Wire **CE**₂₇₉₇

The PTCA guide wire for any case and any fingers

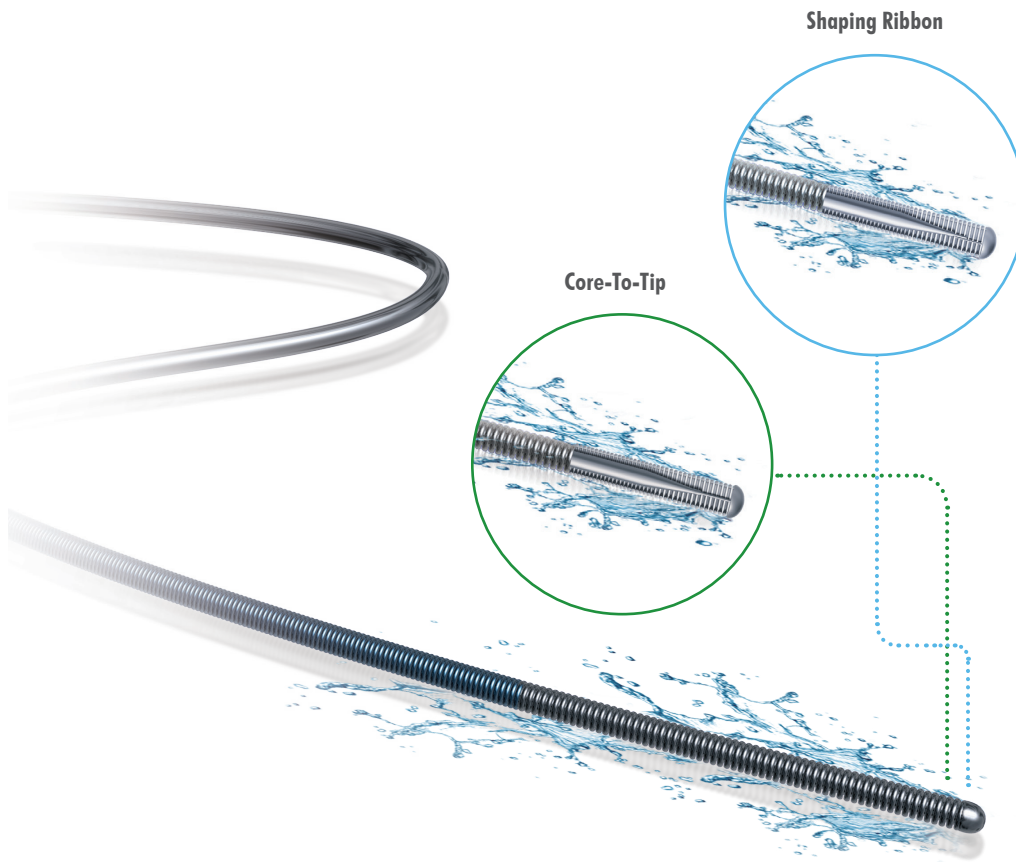

Product Information

Features


- Core-to-tip or Shaping ribbon
- Nitinol + stainless steel core
- Patented welding technique
- Hydroskin™ hydrophilic coating
- 1:1 transmission
- 3 tip-load levels

Benefits

- Multiple choices
- Good visibility
- Super-smooth navigation
- Excellent tactile feedback
- Balanced flexibility and support
- Great pushability and trackability

Ranges

- O.D: 0.014"
- Length: 185/300cm
- Tip load: floppy/supersoft/soft
- Tip shape: straight/J-tip
- Core-to-tip / shaping ribbon
- Radial support: moderate

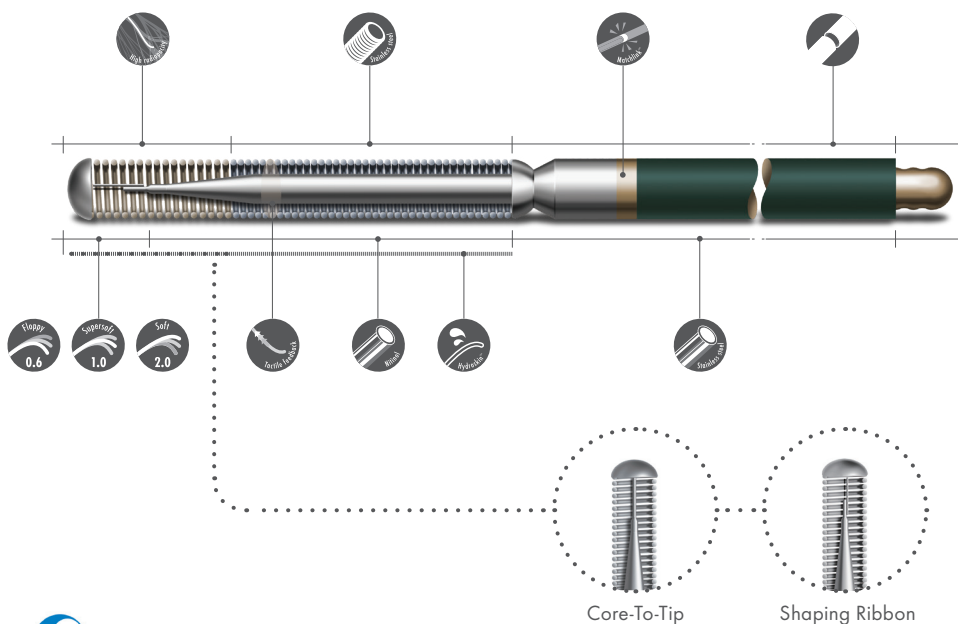

- 3 Tip-load Levels offer multiple choices for you
- 3 cm Platinum Coil provides great visibility
- Stainless Steel Spring Coil balances flexibility and support
- ACT™ Welding enhances tactile feedback
- 30 cm Hydroskin™ Coating guarantees smooth navigation
- Nitinol Core Wire improves flexibility/trackability
- Matchlink™ Welding ensures 1:1 transmission
- Stainless Steel Core Wire provides good pushability
- PTFE Coating reduces friction in catheter

Ordering Information

Core-To-Tip

Catalogue No.	O.D.	Length	Tip Shape	Tip Load	Radial Support
GW B185-0	0.014"	185cm	Straight	Floppy (0.6g)	Moderate
GW B185-0J	0.014"	185cm	J	Floppy (0.6g)	Moderate
GW B300-0	0.014"	300cm	Straight	Floppy (0.6g)	Moderate
GW B300-0J	0.014"	300cm	J	Floppy (0.6g)	Moderate
GW B185-1	0.014"	185cm	Straight	Supersoft (1g)	Moderate
GW B185-1J	0.014"	185cm	J	Supersoft (1g)	Moderate
GW B300-1	0.014"	300cm	Straight	Supersoft (1g)	Moderate
GW B300-1J	0.014"	300cm	J	Supersoft (1g)	Moderate
GW B185-2	0.014"	185cm	Straight	Soft (2g)	Moderate
GW B185-2J	0.014"	185cm	J	Soft (2g)	Moderate
GW B300-2	0.014"	300cm	Straight	Soft (2g)	Moderate
GW B300-2J	0.014"	300cm	J	Soft (2g)	Moderate

Shaping Ribbon

Catalogue No.	O.D.	Length	Tip Shape	Tip Load	Radial Support
GW B185-S0	0.014"	185cm	Straight	Floppy (0.6g)	Moderate
GW B185-S0J	0.014"	185cm	J	Floppy (0.6g)	Moderate
GW B300-S0	0.014"	300cm	Straight	Floppy (0.6g)	Moderate
GW B300-S0J	0.014"	300cm	J	Floppy (0.6g)	Moderate
GW B185-S1	0.014"	185cm	Straight	Supersoft (1g)	Moderate
GW B185-S1J	0.014"	185cm	J	Supersoft (1g)	Moderate
GW B300-S1	0.014"	300cm	Straight	Supersoft (1g)	Moderate
GW B300-S1J	0.014"	300cm	J	Supersoft (1g)	Moderate
GW B185-S2	0.014"	185cm	Straight	Soft (2g)	Moderate
GW B185-S2J	0.014"	185cm	J	Soft (2g)	Moderate
GW B300-S2	0.014"	300cm	Straight	Soft (2g)	Moderate
GW B300-S2J	0.014"	300cm	J	Soft (2g)	Moderate

Core-To-Tip

Shaping Ribbon

